



AB GERMA®

ALL IN ONE VACUUM
MATTRESS



PRODUCT SHEET

The Germa All in One vacuum mattress is our original high quality, long-life model.

The vacuum mattress is a combined stretcher and vacuum mattress. To save storage space the mattress can lie on the stretcher at all times, acting as a stretcher mattress. If vacuum immobilisation is needed, the vacuum mattress is integrated in the stretcher mattress and is ready for use immediately. Of course, the mattress can also be used as a stand-alone ordinary vacuum mattress.

PRODUCT SPECIFICATIONS

Approved according to European Standard
EN 1865 and EN 1789
MDR 2017/745 Class I Medical device

PRODUCT INFO

- 6 handles on each side – carry capacity 250 Kg
- One-piece back, providing the patient with even spine support.
- unique chambers system

Warranty period: 3 years

MATERIAL

- Outer material: Polyester fibre, PVC-coated, flame retardant
- Inner bag: Polyester
- Granules: Polystyrene granules (flame retardant)


SPECIFICATIONS

- Length: 210 cm
- Width: 115 cm
- Height: 5,5 cm
- Weight 8,7 kg with restraint system
- Loading capacity: 250 kg

TEST SPECIFICATIONS

- Maximum shrinkage on length: 1%
- Maximum shrinkage on width: 3%
- Working temperature of the mattress is -30°C to +70°C

ALL IN ONE VACUUM MATTRESS

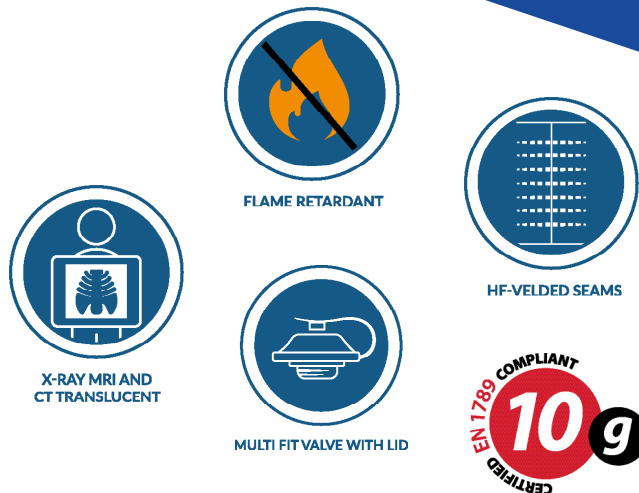

ORDER INFORMATION

Art no: 271101000
Vacuum Mattress AIO 115CM

Art no: 271101001
Vacuum Mattress AIO 115 cm Ready2Go Kit
– including restraint system, footrest Vacuum Pump, repair kit, bag
Weight: 10,6 kg

Art no: 271102000
Vacuum Mattress AIO 100CM

Art no: 271102001
Vacuum Mattress AIO 100CM Ready2Go Kit

Art no: 271103000
Vacuum Mattress AIO 90CM

Art no: 271103001
Vacuum Mattress AIO 90CM Ready2Go Kit

SPARE PARTS

Art no: 12633800000
Art: Germa vacuum valve Black

Art no: 271000501
Art: Repair Kit Mattresses/Splints

Visiting address
AB Germa
Industrig. 54-56
SE-291 36 Kristianstad, Sweden

Post address
AB Germa
Box 2127
SE-291 02 Kristianstad, Sweden

www.germa.se
Phone: +46 (0)44 12 30 30
info@germa.se

BENEFITS

Made to correspondent to EN 1789 and 1865
Approved for European standard and crash test for 10-G

One-piece back, providing the patient with even spine support. Also allows preliminary X-ray diagnostics, direct in to x-ray no extra handling of patient.

Gives rescuers quick handling and excellent immobilisation for the patient.

Deliver consistent and rigid support to the patient safe and convenient immobilisation.

Handles feature an ergonomic grip, making patient lifts as comfortable as possible for rescuers. Handles are reinforced so a patient can be lifted in a patient lifting device.

Handles fixed to the bottom of mattress, maintain immobilisation when lifting patient.

Robust welding—fully welded without any open seams or textiles, making hygienic cleaning of the mattress easy.

Valve with special designed hole, ensuring almost all suction devices fit directly into the valve, no adapters are needed. Reverse air lock system, no closing of valve at stressful site. Made from plastic makes it ideal for X-ray.

Improved unique granules-filled chambered inner bag, ensuring the granules is evenly dispersed to give the patient stable and optimal immobilisation.

Safety belt fixation close to the mattress, safe and comfortable for patient.

PVC coated fabric makes it easy to clean and disinfect with most common chemicals.

Outer material is made of flame retardant PVC-coated fabric. It is flexible making it easy to mould the mattress around the patient.

ALL IN ONE VACUUM MATTRESS

Robust and efficient material to hold vacuum for more than 24 hours, last longer and safer for patient make trustfully immobilisation.

High load capacity 250 kg, safe to handle bariatric patients.

Two different closing mechanisms on the straps, either Hooks & Loops or plastic buckle, allowing the user to choose according to preference.

Repair Kit easy to use for effective repair, no return needed.

Folds easily, in case the stretcher mattress facility is not required.

High quality footrest aluminium vacuum pump evacuate with more than – 700 hpa, quick and time saving.