




AB GERMA®

HEAVY DUTY TRANSFER MATTRESS

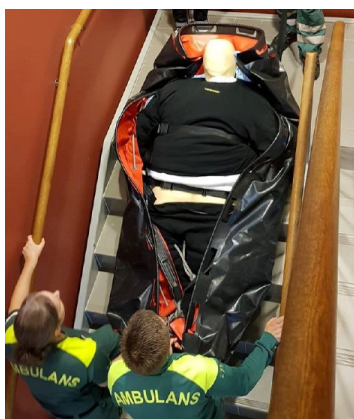

PRODUCT SHEET

Heavy Duty Transfer Mattress is developed to facilitate heavy duty patient handling. The mattress consist of durable and strong PVC coated polyester fabric material equipped with 6 carry handles on each long side plus evacuation handles on head and foot end of the mattress. This equipment is specifically designed to support and help the paramedics and rescuers to be able to transfer and transport a heavy duty patient. The construction and design make it possible to slide the mattress on the floor or down stair cases with patient onboard the mattress. When needed several rescuers are able to lift and carry the patient over obstacles or on to a heavy duty stretcher by using the handles. In next step, from the stretcher to a hospital bed. The purpose is to support handling of heavy duty patient with safety and improved ergonomics.

HEAVY DUTY TRANSFER MATTRESS



PRODUCT SPECIFICATIONS

MDR 2017/745 Class 1 Medical device

PRODUCT INFO

Warranty period: 2 year

- 6 carry handles on each side make it easy to find balance for lifting when needed and carry a patient on board. Easy accessible handles up from ground lift.
- 2 handles on head and foot end to slide the patient on floors and staircases
- 2 Safety belts to secure the patient and cover sheet to protect and hold the patient
- Durable and strong material and design with a carry capacity of 380 kg.

MATERIAL

Material: PVC—polyester fabric material (flame retardant)

PP-plate and grass fibre material to reduce friction

Rope: Nylon

Foam: 60 mm, Polyurethane foam

SPECIFICATIONS SIZE:

Length: 205 cm

Width: 130 cm


HEIGHT: 60 mm

WEIGHT: 12,8 kg

ORDER INFORMATION:

Art: Heavy-Duty Tranfer Mattress

Art no: 271050000

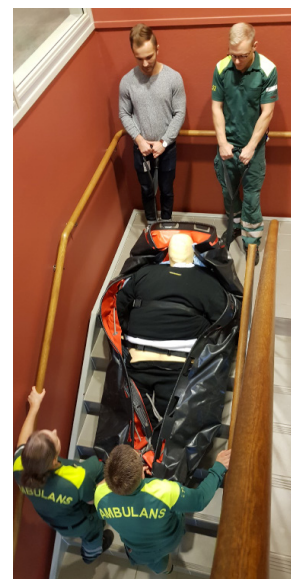

Visiting address
AB Germa
Industrig. 54-56
SE-291 36 Kristianstad, Sweden

Postal address
AB Germa
Box 2127
SE-291 02 Kristianstad, Sweden

www.germa.se
Phone: +46 (0)44 12 30 30
info@germa.se

BENEFITS

MDR 2017/745 Class 1 Medical device

One-piece back, providing the patient with even spine support

Also allows preliminary X-ray diagnostics, direct in to x-ray no extra handling of patient.

Ergonomic handles, making patient lifts as comfortable as possible for rescuers

Handles are reinforced so a patient can be lifted in another patient lifting device.

Handles fixed to the bottom of mattress maintain stable when lifting the patient.

Robust welding—fully welded without any open seams or textiles, making hygienic cleaning of the mattress easy.

Outer material is made of flame retardant PVC-coated fabric.

Extra layer built under the mattress, safety, minimize risk of damage the fixation, last longer.

High load capacity 380 kg, safe to handle heavy duty people.

Closing mechanisms of the straps, plastic buckle Fix lock, allowing the user to quick tighten and adjust restraints on the patient.

HEAVY DUTY TRANSFER MATTRESS

PVC coated fabric makes it easy to clean and disinfect with most common disinfectants, easy to keep clean and ready for use.

Repair Kit easy to use for effective repair, quick to repair no return needed.

Visiting address
AB Germa
Industrig. 54-56
SE-291 36 Kristianstad, Sweden

Postal address
AB Germa
Box 2127
SE-291 02 Kristianstad, Sweden

www.germa.se
Phone: +46 (0)44 12 30 30
info@germa.se