



AB GERMA®

RESCUE LITE FULL BODY SPLINT



PRODUCT SHEET

Intended use is to give good stabilization for the patient being rescued in combination with e.g. a basket stretcher. Light construction for rescuers to carry with them in difficult terrain.

Rescue Lite Full Body Splint design is based on feedback from medical and rescue experts with great experience from rescue operations in difficult terrain. Construction based on original high quality materials traditionally used producing Germa world class vacuum mattresses.

Rescue Lite Full Body Splint is utilizing improved body shape design to optimize the vacuum full body splint capabilities. Offers 35% saving on storage space compared to Germa Easyfix Plus II Vacuum Mattress.

RESCUE LITE FULL BODY SPLINT

The Rescue Lite Full Body Splint features a unique design that provides effective stabilization of the entire body, specifically the upper back and neck, while excluding the cervical collar. It also includes a reinforced base plate that covers the entire cervical spine. This together with fixation restraint system using X-system to support and create the most effective immobilisation.




FLAME RETARDANT



X-RAY MRI AND
CT TRANSLUCENT



MULTI FIT VALVE WITH LID

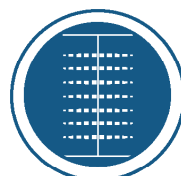

HF-VELDED SEAMS



REINFORCED WOOD
INNER SYSTEM

PRODUCT SPECIFICATIONS

MDR 2017/745 CLASS 1 MEDICAL DEVICE

PRODUCT INFO

- 3 handles on each side
- Removable/changable glass fiber back side
- Inner bag with chamber system
- Carry capacity: 220 kg

Warranty period: 3 years

MATERIAL

- Outer material: Polyester fibre, PVC-coated, flame retardant.
- Inner bag: Polyester
- Granules: Polystyrene granules (flame retardant)
- Head/neck part: Birch plywood

SPECIFICATION

- Length: 200 cm
- Width: 136 (Widest part)/74 (Head part) cm
- Height: 3,5 cm
- Weight: 8,5 kg
- Weight R2G kit: 10 kg

TEST SPECIFICATIONS

- Working temperature of the mattress is -30°C to +70°C

ORDER INFORMATION

Art no: 271513001

RESCUE LITE FULL BODY VACUUM SPLINT

READY2GO KIT, footrest Vacuum Pump, repair kit, Back Pack bag.

SPARE PARTS

Art no: 12633800000

Art: Germa vacuum valve Black

Art no: 15000000112

Art: Repair Ki

Visiting address
AB Germa
Industrig. 54-56
SE-291 36 Kristianstad, Sweden

Post address
AB Germa
Box 2127
SE-291 02 Kristianstad, Sweden

www.germa.se
Phone: +46 (0)44 12 30 30
info@germa.se

BENEFITS

Allows preliminary X-ray & MRI diagnostics, direct in to x-ray & MRI machines with no extra handling of patient.

Glass fiber back side to improve durability and ability to slide the mattress on the ground/floor/stairs etc.

Head, neck and upper back part with wood part to secure stronger immobilisation.

Germa improved chamber technology for quick and easy preparation. Give rescuers quick handling and excellent immobilisation for the patient ensuring the granulate is evenly dispersed to give the patient stable and optimal immobilisation.

Easy to lift patient in/out from Basket Stretcher

Light weighted and small storage, folds easily, 30% less storage space required compared to standard vacuum mattress. Easy to carry in the specially designed back pack.

Robust welding, fully welded without any open seams or textiles. Making hygienic cleaning of the mattress easy.

Valve with special designed hole, ensuring almost all suction devices fit directly into the valve, no adapters are needed. Reverse air lock system, no closing of valve at stressful site.

V-shaped design to make sure upper wide part will make all sized patients the very best immobilisation of the back, neck and head.

X-designed restraint belt fixation close to the patient, support immobilisation process and comfortable for patient. Supplied with five straps with fix lock, in different colours, helping the rescuer to place the straps correctly and easy to position and adjust.

Robust and efficient material to hold vacuum for more than 24 hours, safety for patient and make confiding immobilisation.

Load capacity 220 kg, safe to handle heavy duty people.

RESCUE LITE FULL BODY SPLINT

PVC coated fabric make it easy to clean and disinfect with most common disinfectants.

Repair Kit easy to use for effective repair, quick to repair no return needed.

Deliver consistent and good immobilisation of the patient in combination with a Basket stretcher.

High quality footrest aluminium vacuum pump evacuate with more than -700 hpa, durable, quick and time saving.