




AB GERMA®

TRANSFER MATTRESS NATO MOD STRETCHER

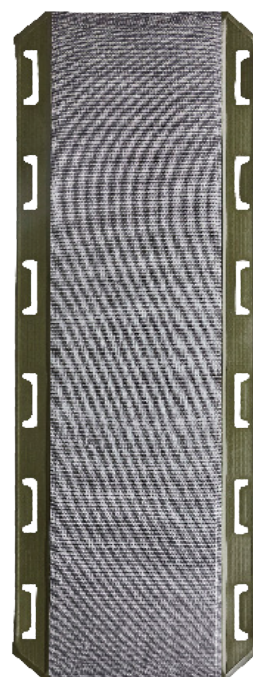

PRODUCT SHEET

Military paramedics needed a Transfer Mattress MOD that could not only be used as a slidable transfer mattress as well as a comfort mattress on a NATO stretcher but could also be used in conjunction with the NATO stretcher belts system for secure transports. We modified our popular standard Transfer Mattress 12 Handles MOD with slots where the belts can be pulled through for secure transport on the NATO Stretcher.

In all other aspects it is as durable, strong, slidable and comfortable as our standard Transfer Mattress 12 Handles MOD.

TRANSFER MATTRESS NATO MOD STRETCHER



X-RAY MRI AND
CT TRANSLUCENT



FLAME RETARDANT



HF-VELDED SEAMS

PRODUCT SPECIFICATIONS

MDR 2017/745 Class 1 Medical device

PRODUCT INFO

- 6 handles on each side
- Glass fiber back side
- Slots in the mattress for the stretcher belt system

Warranty period: 3 years

MATERIAL

- Outer material: Polyester fibre, PVC-coated, flame retardant and Non IR reflective
- Inner material: 30 mm thick Polyurethane Foam

SPECIFICATION

- Length: 187 cm
- Width: 65 cm
- Height: 3,0 cm
- Weight 4,5 kg

ORDER INFORMATION

Art no: 271008881

Transfer Mattress MOD NATO Stretcher

SPARE PARTS

Art no: 15005002817

Art: Repair Kit

Visiting address
AB Germa
Industrig. 54-56
SE-291 36 Kristianstad, Sweden

Post address
AB Germa
Box 2127
SE-291 02 Kristianstad, Sweden

www.germa.se
Phone: +46 (0)44 12 30 30
info@germa.se

BENEFITS

MDR 2017/745 Class 1 Medical device

TRANSFER MATTRESS NATO MOD STRETCHER

Easy to carry or slide over a body on the mattress with 12 strong handles, 6 on each long side.

Handle parts can be folded under the mattress on the stretcher.

Sliding surface on the reverse of fiberglass material, reducing friction and allow sliding of the mattress e.g. on the floor/ground or from stretcher to stretcher.

Should be used with restraints to secure the patient for safety reasons.

Comfort, ethical and ergonomic handling of the patient.

Hygienic; fully welded and easy to clean with most common disinfectants.

Fits together with Germa Step by Step evacuation Kit